

Review of Vatech SW Integration with AXS Platform

vatech

SW & Service Overview

EzDent-i

- 2D Viewer: Select patient, initiates scan session via IO Scanner Scan SW, and manages imaging data

Clever One

- Integrated 2D & 3D Viewer : Select patient, initiates scan session via IO Scanner Scan SW, and manages imaging data

EzServer

- On-premise central DB Server : Data Storage & Management

Clever Space

- Cloud Service for Imaging Data Sharing

Clever Lab

- Service for Dental Lab Operations Management

Data Integration Prerequisites

Data Management Overview

- All data within the dental clinic is stored in EzServer (intraoral scanner integration is mandatory [**Under discussion*])
- Dental lab orders are transmitted to Clever Lab


Unidirectional Data Flow

- Only one-way transmission in the direction of EzServer → AXS → Clever Lab is supported
- Bidirectional synchronization (Sync) between AXS ↔ EzServer is not supported

Data Integrity and History Preservation

- Once a case is saved, it cannot be overwritten; a new Case ID must be issued for any modifications
- Linkage between the original case and the modified case must be maintained (history tracing must be supported)

Basic Workflow (EzServer → AXS → Clever Lab)

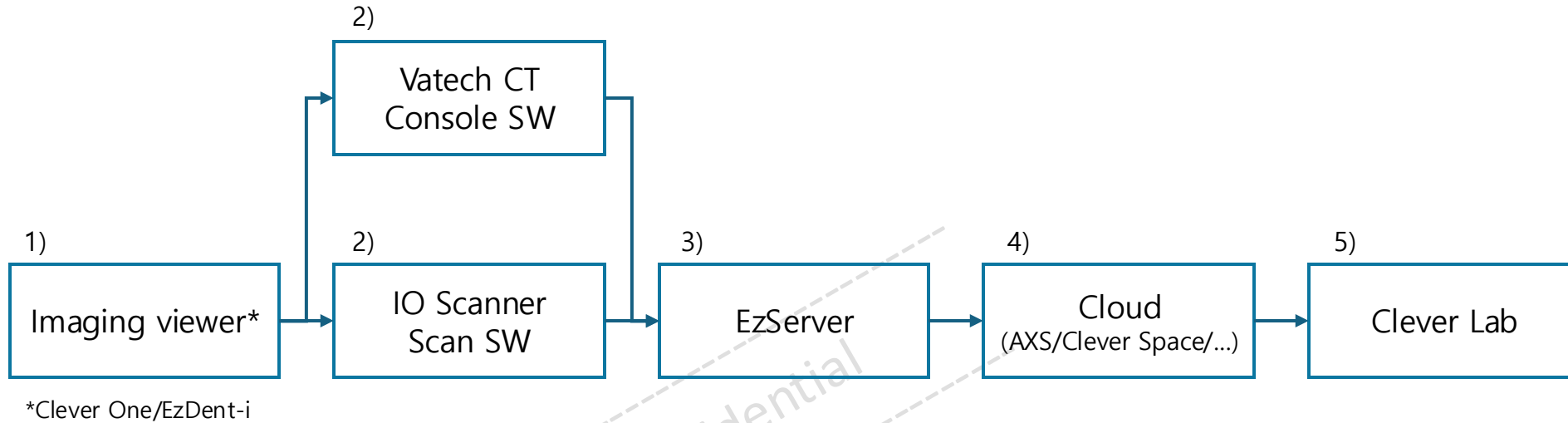


① → ② : Launch Scan SW and Console SW

② → ③ : Save Scan and CT files

③ → ④ : Upload patient information, Image Data, Case Information

④ → ⑤ : Transfer patient information, Image Data, Case Information


Integration Scenario 1. Basic Data Flow (EzServer → AXS)

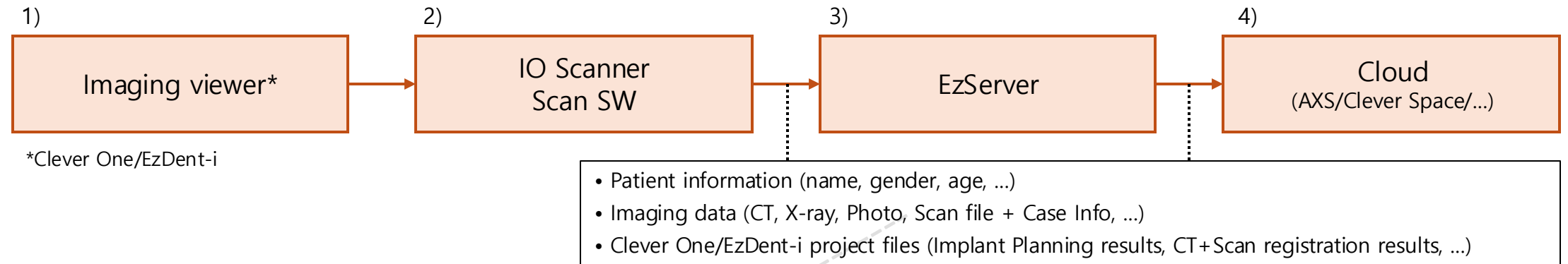

Step	Stage	Description	Notes
1	Initiate scan from Imaging Viewer	Select a patient in the Imaging Viewer and initiate an intraoral scan → IO Scanner Scan SW creates a new case based on registered patient information and prepares for scanning based on Case ID	
2	Scan and generate/save scan file	Save original scan file + Case Info locally	
3	Save scan file to EzServer	Automatically transmit locally saved scan file (from Step 2) from Imaging viewer to EzServer	Integration method to be discussed
4	Upload to AXS	One-way transmission: EzServer → AXS	

Integration Scenario 2. Modifying Data (Saved as New Case)


(Saved as New Case When Rescanning Due to Scan Quality Issues, etc.)

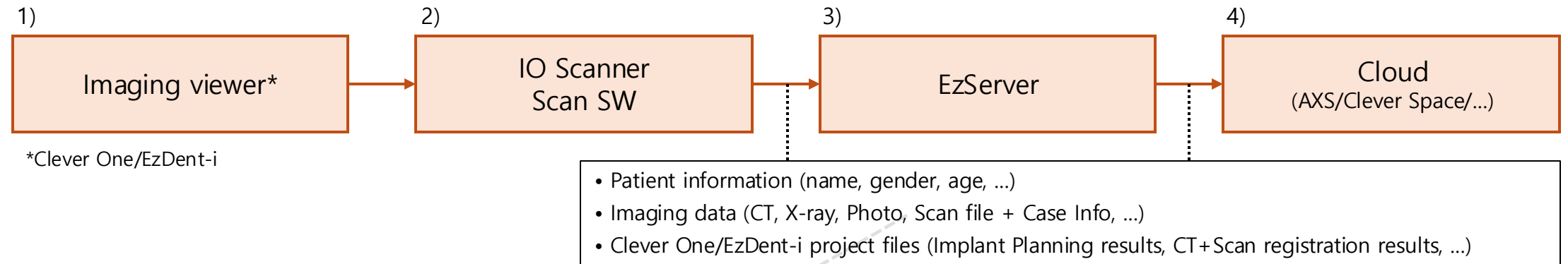

Step	Stage	Description	Notes
1	Select case to modify	Select patient in Imaging viewer → pass chart number + Case ID to IO Scanner Scan SW	
2	Modify existing scan data	Launch Scan SW with patient and case pre-selected based on chart number + Case ID → user modifies scan data (rescan or edit) → save locally under new Case ID	Modified data saved separately under new Case ID without overwriting existing case
3	Save scan file to EzServer	Automatically transmit locally saved scan file (from Step 2) from Imaging viewer to EzServer	
4	Upload to Cloud	One-way transmission: EzServer → Cloud	

Integration Scenario 3. AXS - Clever Lab

(Scenario 1: Order Information Transmission.)

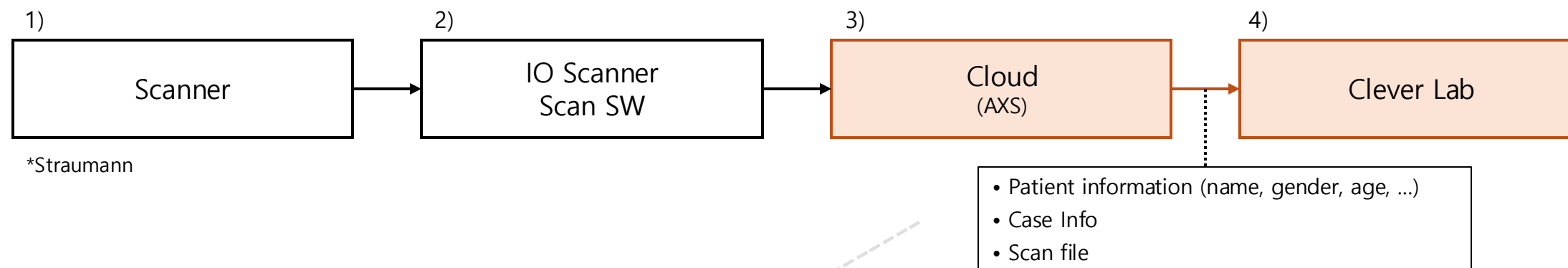

Step	Stage	Description	Notes
1	Patient Selection & Order Creation	Select patient in Straumann Scan SW → Enter prosthetic information (tooth number, prosthesis type, material, etc.) and create order	Unique Order ID assigned upon order creation
2	Scan & File Generation	Perform intraoral scan with Straumann Scanner → Scan file + order information saved in Straumann Scan SW(Server)	Case ID and Order ID mapped
3	Data Transmission to AXS	Straumann Scan SW → AXS: transmit scan file + order information	Transmission method (API/file transfer, etc.) to be discussed
4	Clever Lab Order Registration	AXS → Clever Lab: deliver order information → Order automatically registered in Clever Lab and status set to pending	Data mapped based on Order ID

Integration Scenario 3. AXS - Clever Lab

(Scenario 2: Work Status Synchronization.)



Step	Stage	Description	Notes
1	Work Status Update	Update work status in Clever Lab (Order Received → In Progress → Completed, etc.)	Timestamp recorded upon status change
2	Status Transmission to AXS	Clever Lab → AXS: transmit updated work status	Status synchronized based on Order ID

Integration Scenario 3. AXS - Clever Lab

(Scenario 3: Confirmation Request)



Step	Stage	Description	Notes
1	Confirmation Request Submission	Clever Lab completes work → Creates confirmation request with design file or try-in result attached	Unique Confirm ID assigned upon confirmation request
2	Delivery via AXS	Clever Lab → AXS → Straumann Console: deliver confirmation request data	Design images, scan files, etc. can be attached
3	Confirmation Result Transmission	AXS → Clever Lab: deliver confirmation result	Result mapped based on Confirm ID
4	Clever Lab Follow-up	If approved: proceed to fabrication / If revision requested: rework and resubmit confirmation request	Scenario 3 can be repeated

End of Document

Thank you

vatech